

QuDAP: an Open Source Software for Quantum Materials Characterization

Skai White¹, Jingyu Jia², Chunli Tang³, Masoud Mahjouri-Samani³, and Wencan Jin³

¹ Department of Mathematics, Hampton University, Hampton, VA, USA ² Great Neck South High School, Lake Success, NY, USA ³ Auburn University, Auburn, AL, USA

DOI: [10.xxxxxx/draft](https://doi.org/10.xxxxxx/draft)

Software

- [Review](#)
- [Repository](#)
- [Archive](#)

Editor: [Open Journals](#)

Reviewers:

- [@openjournals](#)

Submitted: 01 January 1970

Published: unpublished

License

Authors of papers retain copyright and release the work under a Creative Commons Attribution 4.0 International License ([CC BY 4.0](#)).

Summary

Quantum Materials Data Acquisition and Processing (QuDAP), a Python-based and open-source software package, is designed to control and automate material characterizations based on the Physical Property Measurement System (PPMS) ([Quantum Design, n.d.](#)). The software supports major hardware interfaces and protocols (USB, RS232, GPIB, and Ethernet), enabling communication with the measurement modules associated with the PPMS. It integrates multiple Python libraries to realize instrument control, data acquisition, and real-time data visualization. Here, we present features of QuDAP, including direct control of instruments without relying on proprietary software, real-time data plotting for immediate verification and analysis, full automation of data acquisition and storage, and real-time notifications of experiment status and errors. These capabilities enhance experimental efficiency, reliability, and reproducibility.

Statement of need

The exploration of quantum materials ([Awschalom et al., 2021](#); [Azam et al., 2022](#); [Kumari et al., 2021](#); [Wei et al., 2020](#)) has gained significant interest over decades, benefiting from their unique properties at the quantum level, compact size, and low energy consumption. These materials are promising candidates for applications in future quantum computers, sensors, and memory devices. Studying the intriguing properties of quantum materials demands extreme conditions such as low temperature, high magnetic field, and ultra-high vacuum environments. Therefore, many experimental approaches have been employed to characterize quantum materials, including magnetometry ([Foner, 1956](#)), spectroscopic techniques ([Rahman et al., 2021](#); [Tang et al., 2023](#)), electrical transport ([Nagaosa et al., 2010](#); [Rojo et al., 2013](#)), thermal transport ([Goyal et al., 2023](#); [Kalantari & Zhang, 2022](#)), and advanced microscopy ([Li et al., 2018](#); [Zhang et al., 2018](#)). Among them, the Physical Property Measurement System (PPMS, Quantum Design) provides the capability to achieve those goals with precise control of temperature and magnetic field under a high vacuum environment. In this process, PPMS needs to communicate with multiple peripheral instruments through hardware interfaces, including USB, RS232, GPIB, and Ethernet. To control and automate the data acquisition, software is needed to integrate the Standard Commands for Programmable Instruments (SCPI) protocol ([Keysight Technologies, 2022](#)).

QuDAP is a Python-based data acquisition and processing software package to control and automate quantum material characterization experiments. It is developed based on the measurement modules of Quantum Design DynaCool PPMS including magnetometry, electrical transport, and ferromagnetic resonance spectroscopy. This software package provides a graphic user interface to visualize and process the experimental data in real-time, giving the users a

more efficient way to monitor the experiment. [Figure 1](#) demonstrates the connection utility of this software. All the instruments are connected through the hardware interface of USB, RS232, GPIB, and ethernet. The commands for each instrument are sent using PyVISA ([Grecco et al., 2023](#)) and MultiPyVu ([Damon D Jackson, 2023](#)). The software provides the benefits as summarized below: 1. Provide direct Python script communication and control of PPMS and instruments without using the built-in software, which improves the tunability and efficiency of the experiment. 2. Built-in demagnetization process before each measurement to enhance the reliability of the measurement. 3. Fully automated data acquisition and saving process with real-time plotting and progress visualization. 4. Save the data with specific identifiers to avoid data overwrite and record the experiment configuration of each measurement. Real-time notification on the measurement status and program error through push notification, allowing the user to promptly identify and verify the experimental and parameter setup. 5. Various post data processing features allowing to process optical, magnetic, and electrical data.

In the following sections, we introduce the peripheral instruments we have integrated into our software and show the experiment procedure for this software.

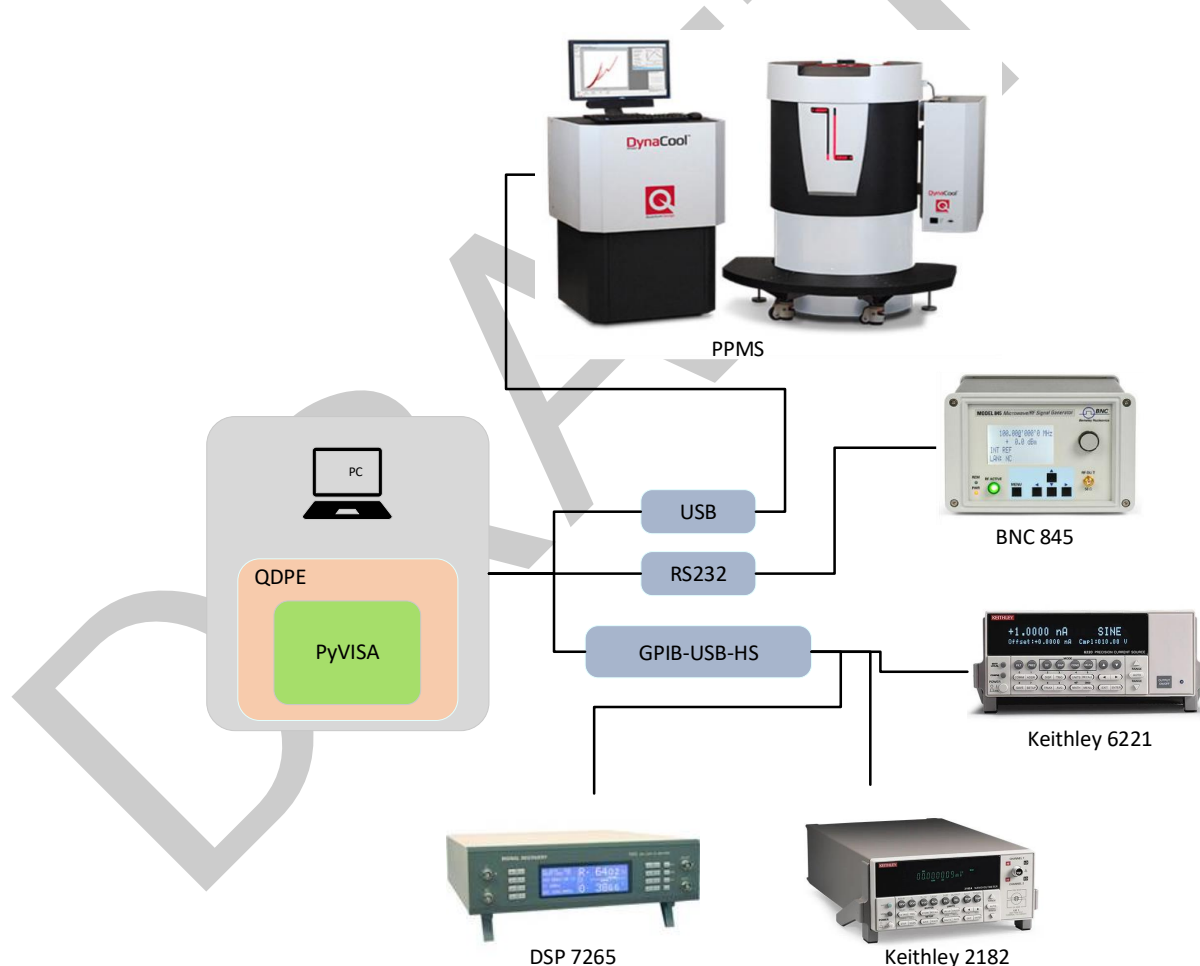


Figure 1: Block Diagram of the instrumentation of QuDAP. The QuDAP software establish the instrument communication using USB, RS232, and GPIB to control the PPMS ([Quantum Design, n.d.](#)), BNC845 ([Berkeley Nucleonics, n.d.](#)), Keithley 6221 ([Tektronix, n.d.-a](#)), Keithley 2182 ([Tektronix, n.d.-b](#)), and DSP7265 ([Signal Recovery, n.d.](#)).

57 Graphic User Interface

58 QuDAP utilizes MultiPyVu for controlling the PPMS. It is a socket-based server/client that
59 enables both readback and control of the temperature, magnetic field, and chamber operations
60 of the sample environment. MultiPyVu also contains a wait function that allows the user
61 to pause the operation until the desired stability criteria are met, while an additional class
62 simplifies the writing of MultiVu-readable data files.

63 The interface of the software is shown in [Figure 2](#). It contains three parts from left to right:
64 main menu, options sidebar, and run window. The main menu allows users to select the
65 option of the dashboard, data processing, and experiment. The detailed options are shown
66 in the options sidebar once select the main option. The software starts up to the dashboard
67 window, which shows all the connections from VISA. The dashboard window also provides the
68 feature of quick access to the functions and directs to the selected page. The data processing
69 options feature data processing from various quantum materials characterization, including
70 magnetometry, ferromagnetic resonance, and electrical transport. The last option on the main
71 menu is the experiment option. Under the experiment option, there are features of individual
72 control for each instrument, preset experiment type, and quick test. The individual instrument
73 gives the user many options to control each feature. Besides, QuDAP allows users to customize
74 the experiment setup to link and control instruments. The preset measurements can be found
75 under the Measurement tab in the option sidebar. This feature provides the users with the
76 preset script for both ETO and FMR measurements. Once selecting the preset option, the user
77 will connect instruments and set up the parameters. A log panel will pop up when the user
78 presses the start button and the experimental log will be recorded as a txt file. To facilitate the
79 data management, an example of the data file naming convention will be provided. After the
80 measurement starts, in the backend, the software will check all the parameters to ensure a valid
81 entry. The real-time data point will be shown in the run window and users will receive push
82 notifications from Telegram once the measurement is done or some unexpected error happens.
83 The quick Test option allows the users to do the measurements without the constraint of the
84 preset script. Users can select the instruments they need for specific experiments. After the
85 measurements are finished, users can find metadata as well as the recorded experimental under
86 the preselected folder and process the data under the data processing option.

87 This software is available in the format of a Python script file, and users can run the program
88 based on their operating system. It will also pack into an executable software for those who
89 do not have experience in Python programs.

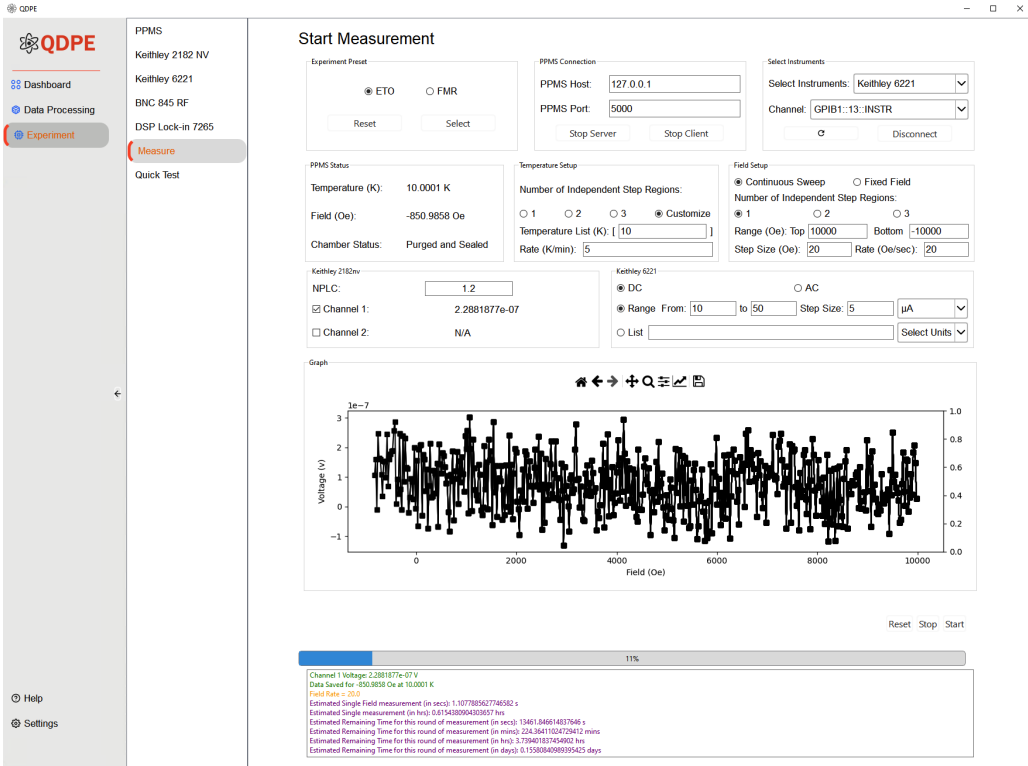


Figure 2: GUI Interface of the QuDAP. The interface shows the ETO measurement setup and real-time data visualization.

Acknowledgements

S. W was funded by Collaborative Approaches among Scientists and Engineers REU National Science Foundation EEC Division Grant 2349639.

References

Awschalom, D. D., Du, C. R., He, R., Heremans, F. J., Hoffmann, A., Hou, J., Kurebayashi, H., Li, Y., Liu, L., Novosad, V., & others. (2021). Quantum engineering with hybrid magnonic systems and materials. *IEEE Transactions on Quantum Engineering*, 2, 1–36.

Azam, N., Boebinger, M. G., Jaiswal, S., Unocic, R. R., Fathi-Hafshejani, P., & Mahjouri-Samani, M. (2022). Laser-assisted synthesis of monolayer 2D MoSe2 crystals with tunable vacancy concentrations: Implications for gas and biosensing. *ACS Applied Nano Materials*, 5(7), 9129–9139.

Berkeley Nucleonics. (n.d.). *Model 845 datasheet v2.75*. https://www.berkeleyelectronics.com/sites/default/files/products/resources/845_datasheet_v275.1.pdf

Damon D Jackson, Q. D. (2023). *MultiPyVu 2.1.4*. <https://pypi.org/project/MultiPyVu/>

Foner, S. (1956). Vibrating sample magnetometer. *Review of Scientific Instruments*, 27(7), 548–548.

Goyal, H., Peek, S. E., Sellers, J. A., & Hamilton, M. C. (2023). Methodology to characterize thermal properties of thin film superconductors using a DynaCool physical property measurement system. *IEEE Transactions on Applied Superconductivity*, 33(5), 1–5.

- 109 Grecco, H. E., Dartailh, M. C., Thalhammer-Thurner, G., Bronger, T., & Bauer, F. (2023).
 110 PyVISA: The python instrumentation package. *Journal of Open Source Software*, 8(84),
 111 5304. <https://doi.org/10.21105/joss.05304>
- 112 Kalantari, M. H., & Zhang, X. (2022). Thermal transport in 2D materials. *Nanomaterials*,
 113 13(1), 117.
- 114 Keysight Technologies. (2022). *Instrument automation with python*.
- 115 Kumari, S., Pradhan, D. K., Pradhan, N. R., & Rack, P. D. (2021). Recent developments on
 116 2D magnetic materials: Challenges and opportunities. *Emergent Materials*, 4(4), 827–846.
- 117 Li, X., Sun, M., Shan, C., Chen, Q., & Wei, X. (2018). Mechanical properties of 2D materials
 118 studied by in situ microscopy techniques. *Advanced Materials Interfaces*, 5(5), 1701246.
- 119 Nagaosa, N., Sinova, J., Onoda, S., MacDonald, A. H., & Ong, N. P. (2010). Anomalous hall
 120 effect. *Reviews of Modern Physics*, 82(2), 1539–1592.
- 121 Quantum Design. (n.d.). *Quantum design PPMS DynaCool*. [https://www.qdusa.com/
 122 products/dynacool.html](https://www.qdusa.com/products/dynacool.html)
- 123 Rahman, S., Torres, J. F., Khan, A. R., & Lu, Y. (2021). Recent developments in van der waals
 124 antiferromagnetic 2D materials: Synthesis, characterization, and device implementation.
 125 *ACS Nano*, 15(11), 17175–17213.
- 126 Rojo, M. M., Calero, O. C., Lopeandia, A., Rodriguez-Viejo, J., & Martin-Gonzalez, M. (2013).
 127 Review on measurement techniques of transport properties of nanowires. *Nanoscale*, 5(23),
 128 11526–11544.
- 129 Signal Recovery. (n.d.). *Model 7265 DSP lock-in amplifier*. [https://www.atecorp.com/
 130 atecorp/media/pdfs/data-sheets/7265lockinampv7web.pdf?ext=.pdf](https://www.atecorp.com/atecorp/media/pdfs/data-sheets/7265lockinampv7web.pdf?ext=.pdf)
- 131 Tang, C., Alahmed, L., Mahdi, M., Xiong, Y., Inman, J., McLaughlin, N. J., Zollitsch, C.,
 132 Kim, T. H., Du, C. R., Kurebayashi, H., & others. (2023). Spin dynamics in van der waals
 133 magnetic systems. *Physics Reports*, 1032, 1–36.
- 134 Tektronix. (n.d.-a). *Keithley 6220-6221*. <https://download.tek.com/datasheet/6220-6221.pdf>
- 135 Tektronix. (n.d.-b). *Model 2182 and 2182A nanovoltmeter user's manual*. [https://download.
 136 tek.com/manual/2182A-900-01_May_2017.pdf](https://download.tek.com/manual/2182A-900-01_May_2017.pdf)
- 137 Wei, S., Liao, X., Wang, C., Li, J., Zhang, H., Zeng, Y.-J., Linghu, J., Jin, H., & Wei, Y. (2020).
 138 Emerging intrinsic magnetism in two-dimensional materials: Theory and applications. *2D
 139 Materials*, 8(1), 012005.
- 140 Zhang, H., Huang, J., Wang, Y., Liu, R., Huai, X., Jiang, J., & Anuso, C. (2018). Atomic
 141 force microscopy for two-dimensional materials: A tutorial review. *Optics Communications*,
 142 406, 3–17.